

Hans Henke

Preparative Low Pressure Liquid Chromatography










For further information and updates regarding this book, please use our online service Info Click



- 1. Go to www.vogel-buchverlag.de
- 2. Click on the button Info Click
- 3. Enter your access code and e-mail address

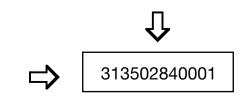

Your access code

Graduate engineer Dr. Hans Henke

Preparative Low Pressure Liquid Chromatography

Vogel Buchverlag

Biographical sketch:

Dr. techn. Dipl.-Ing. Hans Henke, born in 1934, studied technical chemistry at the Technical University of Graz, Austria, from 1956 to 1964. His Diploma Thesis was entitled Quantitative Paper Chromatographic Determination of Vitamin C in Fruit and Vegetables and his Doctoral Thesis addressed the analysis of components of nucleic acids (first quantitative gel Chromatographic determination on Sephadex G gels). During 1964 to 1970, Dr. Henke was Research Assistant at the Institute of Biochemistry, Food Chemistry, and Microchemistry of the Technical University of Graz. From 1970 to 1994 he was group leader of the analytical department (general analysis, gas and liquid chromatography) of the research facility of AKZO Nobel AG at Obernburg/Main, Germany. During this time he authored numerous publications and internal reports on analytical and preparative LC methods.

A brief selection from these publications:

- New column liquid Chromatographic analysis of spinning preparations and fiber extracts
- New analytical methods for polyurethanes
- New and highly accurate preparative and analytical LC separations of alkylpolyglycosides, sugar derivatives, rubber chemicals, etc.

Translated by Anthony J. Rackstraw

Additional information www.vogel-buchverlag.de

ISBN 978-3-8343-3135-9 1st edition 2011 All rights reserved Copyright 2011 by Vogel Business Media GmbH & Co. KG, Würzburg In recognition of the achievements of the brilliant inventor of Column Adsorption Chromatography (1903/1906) Michail Semenovich TSWETT (1872–1919)

Preface

Preparative column liquid chromatography was and remains an important subdiscipline of chromatography, whether as an indispensable component of instrumental analysis in the isolation of unknown compounds for identification or as a gentle separation or purification procedure, or also as a method for the isolation of sensitive natural products.

On mention of preparative column liquid chromatography, the mind immediately turns to the numerous laboratories engaged in industrial research and those in university institutes grappling with many different separation problems, in the gram range and not on a production scale because that is a different world. In particular, these laboratories require straightforward, versatile, and above all relatively inexpensive separation systems. The samples or substance mixtures to be separated can be of synthetic or natural origin. The latter are generally plant extracts.

The predominant role of thin-layer chromatography in qualitative testing in preparative organic laboratories, as well as for plant constituents, has remained unchallenged for many decades. After TLC monitoring of a synthesis, it usually becomes necessary to isolate the desired synthesis products or even the by-products in preparative amounts.

The principal component of all separation units remains the separation system, consisting of the column packing materials and the mobile phase. Preparative column liquid chromatographic set-ups must be equipped with an economical separation system.

A separation system is economical if neat solvents are used as mobile phases and very cheap or long-lived compound packing materials are used. The feasibility of fulfilling both criteria while producing highly selective separation systems will be demonstrated in this book with the aid of numerous examples.

Use of neat solvents as mobile phases has yet another great advantage for the practitioner:

- □ Simple transfer of all analytical TLC separations from a plate to an almost identical preparative system in a column is guaranteed.
- Probably the most important task in preparative column liquid chromatography is to establish the best neat solvent for the sample, i.e. the substance mixture.
- □ The next very important step is to obtain a qualitative impression of the composition of the sample by TLC and/or HPLC.

These qualitative analyses suffice to provide important information about the substance mixture and permit a straightforward and rapid decision about how to solve the separation task. In preparative low-pressure column liquid chromatography the column packing materials are packed in glass columns of various lengths and IDs. The phases used may be pressure-stable silica gels or swellable organic polymers with an average particle size of $\geq 20 \ \mu\text{m}$. The range of available column packing materials is large but only few of them are versatile in use.

This book has been written in keeping with the motto: "By the practitioner for the practitioner"

with the aim of encouraging the user of preparative methods to solve his own separation problems with the aid of simple and readily reproducible examples taken from everyday practice.

A critical remark is in place at this juncture: Preparative separations with unnamed peaks are completely unsuitable as illustrative examples!

This book presents various laboratory-scale preparative separations, which, however, can be readily scaled-up to pilot plant dimensions.

Complex substance mixtures, whether of synthetic or natural origin, can generally be separated according to the following scheme:

- □ Enrichment of the desired products with one of the systems shown and
- □ Separation and isolation of the desired components with the aid of a selective separation system.

The most important attributes of this book are the use of neat solvents as mobile phases for all separations and the unconditional guarantee of transfer of analytical TLC separations to a preparative column.

Separation processes in practice

Liquid chromatography is performed by two different yet complementary procedures: thin-layer and column chromatography. The former is carried out mainly in the analytical mode and the latter both in the analytical and in the preparative mode. The complementary nature of the two procedures follows from the way in which fast analytical TLC separations can serve as an optimization method for preparative column liquid chromatography. Optimization amounts to determining the most suitable separation system for preparative separation of a substance mixture.

Thin-layer chromatography

This is not the place for a detailed account of thin-layer chromatography (TLC) because several specialized monographs are available; instead we shall focus our attention on transfers from a TLC plate to a preparative column.

It is common knowledge that thin-layer chromatographic separations are generally performed with solvent mixtures as developers. In the case of development with solvents containing two or more components, gradients are generated during development on the silica gel layer. This means that the developing solvent contains less and less polar components the higher the front rises.

Thin-layer chromatographic analysis and sample origin, whether essential oils and plant extracts or samples of synthetic origin, are of enormous practical importance in connection with preparative separations.

It is not always easy to transfer analytical TLC separations to preparative columns, but transfer can prove to be straightforward if two conditions are met:

- □ The sample or substance mixture must undergo at least incipient separation on a silica gel plate when developed with a neat solvent.
- □ As a rule it is advisable to pre-sort the sample components by stepwise elution (see Section V.4.), i.e. to perform sample preparation or enrichment.

A TLC separation with dichloromethane on a silica gel plate is generally very helpful because it provides good qualitative information about the polarity distribution of the sample components. Information is often available at an early stage about the qualitative composition of a sample and that is extremely helpful for separation or enrichment of desired products by stepwise elution.

Problem-free, i.e. guaranteed transfer from an analytical silica gel TLC plate to a preparative column is possible on exclusive use of neat developer. Developers of choice are dichloromethane, ethyl acetate, acetone, and methanol, which are also good solvents for numerous substance mixtures of various compositions – an important prerequisite in preparative column-liquid chromatography.

It is understandable that transfer from silica gel plate to silica gel column is still accepted as a matter of course, or not accepted.

Matters become more problematical on transfer of a separation from a silica gel plate to a column packed with the hydroxypropylated dextran gel Sephadex LH-20, because a transfer from an inorganic to an organic gel can hardly be expected to work. In fact, the opposite could repeatedly be proved true in practice, with the enormous additional benefit that Sephadex LH-20 forms separation systems of exceptional selectivity with neat organic solvents.

Erlenbach

Hans Henke

Contents

Pre	eface	•••••			7
I	Intro	oductio	on		21
II	The	oretica	l Part		
	1	Theor	retical pr	inciples of exclusion chromatography	
		1.1	Volum	e ratios in a gel packing	
		1.2	Elution	parameters in exclusion chromatography	
		1.3	Sephad	ex LH-20: Comparison with HPLC and GPC	
	2	Cross	linked de	extran gels	
		2.1	Produc	tion of Sephadex LH-20	
		2.2	Determ	ination of the interstitial volume (V_0)	
			and the	e pore volume (V_i)	
			2.2.1	Methanol	
			2.2.2	Acetone	
			2.2.3	Dichloromethane	
			2.2.4	Ethyl acetate	
			2.2.5	N-Methyl-2-pyrrolidone (NMP)	
			2.2.6	Water	

IV	Colı	ımn Li	quid Chromatographic Separation Systems	
	1		mn packing materials – Phases	
	2		le phases – Eluents	
		2.1	Solvent mixtures	
		2.2	Neat solvents	61
	3	Prepa	rative columns	
	4	Equip	oment	65
		4.1	Pumps	
		4.2	Sample introduction systems	67
		4.3	Detectors	
		4.4	Fraction collectors	
		4.5	Recording devices	71

V	Prac	tical Part 1 – Preparative LC in Practice	73
	1	Sample preparation	73
	2	Analytical – Preparative	74
		, 1	

	2.1		er chromatography – Column chromatography	
	2.2	HPLC – I	Preparative CLC	78
		2.2.1	High-pressure version	78
		2.2.2	Medium- and low-pressure ranges	82
		2.2.3	HPLC – Sephadex LH-20	95
3	Thin	layer chror	natography – Sephadex LH-20	97
	3.1		tures	
	3.2		products	
			Di- and monocarboxylic esters of	
			4,4´-dihydroxybiphenyl	113
			Glycerol monomyristin phenylurethanes	
			Oligomers of 2,2,4-trimethyl-1,2-dihydroquinoline	
	3.3		racts	
			Hop cone CO ₂ extract	
			Green tea	
			Ginkgo biloba leaves	
4	Separ		mplex substance mixtures	
•	4.1		elution without detection – Procedures	
			Hexane – Dichloromethane	
			Dichloromethane – Methanol	
	4.2		preparations for synthetic fibers	
	1.2		Spinning preparations	
			Properties of modern spinning preparations	
			Composition of spinning preparations	
			Column liquid chromatographic analysis of	140
			commercial preparations, components, and	
			fiber extracts	150
			Extraction of polyester and polyamide fibers	
	4.3	Essential	1, 1,	
	4.4	Resins – 1		
5			ods – Extraction processes	
5	5.1		– Synthetic origin	
	5.1 5.2	Samples -	– Synthetic origin – Natural origin	105
	5.2 5.3	1	0	
			n with neat solvents	
1	5.4		nonitoring method	
6	-) – Separation mechanism A	
	6.1		r characteristics – Compound classes	
	6.2	Molecula		
	6.3		nd number of heteroatoms – Functional groups	
			Polarity	
			Positional isomerism	
			Structural isomerism	
		6.3.4	Double bonds	209

		6.3.5 Oxime functions	220
	6.4	Cis/trans isomerism – Geometrical isomerism	223
	6.5	Diastereomerism	225
7	Sepha	dex LH-20 – Separation mechanism B	
	7.1	Centers of asymmetry – Polarity – Pore spectrum	
	7.2	Extremely selective separations	
	7.3	Comparison separations	
8	Sepha	Idex LH-20 – Scale-up	250
9		ıdex LH-20 – Loadability	
10		Idex LH-20 compared to	
	10.1	Normal phase	270
	10.2	Reversed phase	
	10.3	Normal phase – Reversed phase	
	10.4	HP-Cellulofine	
11	Select	ive preparative separation systems	284
	11.1	Sephadex LH-20/Methanol	
	11.2	Sephadex LH-20/Acetone	290
	11.3	Sephadex LH-20/Ethyl acetate	
	11.4	Sephadex LH-20/Dichloromethane	
	11.5	Sephadex LH-20/Water	312
	11.6	Sephadex LH-20/1-Methyl-2-pyrrolidone	324
	11.7	Packing of glass columns with Sephadex LH-20	327
		11.7.1 Preparative separation column assemblies	332
	11.8	Separation systems – Schemes	334
	11.9	Central preparative laboratory	349
12	Exam	ples of separations from practice	
	12.1	Reaction products	351
		12.1.1 Hofmann degradation	351
		12.1.2 Oligomerization of 2,2,4-trimethyl-1,2-	
		dihydroquinoline	353
		12.1.3 Ethoxylation of urea	364
		12.1.4 Antioxidant – Phenol derivatives	370
		12.1.5 2,5-Norbornadiene addition products	374
		12.1.6 Transformation of epoxides into episulfides	380
		12.1.7 Reaction product of carbon disulfide and	
		dimethylamine	382
13	Glyce	rides	384
	13.1	Monoglycerides	384
		13.1.1 Glyceryl monostearate	385
		13.1.2 Analysis of pure glyceryl monostearate	387
	13.2	Monoglycerides - Synthesis and separation	
		13.2.1 Monoglycerides – Diglycerides	392
	13.3	Diglycerides	

Separ	ation of various products	395
14.1	Pentaerythritol tetraesters	395
14.2	Ethoxylated fatty alcohols	398
14.3	Trimethylolpropane-pelargonate mixture	401
Separ	ation of depolymerization products	403
15.1	Polyurethanes	403
15.2	Polycarbonates	407
15.3	Purely aromatic polyesters	412
15.4	Incomplete butanolyzate of polybutanediol adipate (PBA)	422
Enricl	nment – Quantitation	429
16.1	Quantification of traces of phenol in bisphenol A	430
16.2	Diphenyl ether in polyglycolic acid depolymerization	
10.2	Dipitentyl ether in polygiyeone acte depolymenzation	
10.2	products	431
16.3	1 / 1 /0/ 1 /	
16.3	products	433
16.3 Exam	products	433 434
16.3 Exam	products PAH sample preparation ples of transfers – TLC separation to preparative column From TLC to silica gel column	433 434 435
16.3 Exam 17.1	products PAH sample preparation ples of transfers – TLC separation to preparative column From TLC to silica gel column	433 434 435
16.3 Exam 17.1 17.2	products PAH sample preparation ples of transfers – TLC separation to preparative column From TLC to silica gel column From TLC to Sephadex LH-20	433 434 435 436
16.3 Exam 17.1 17.2 17.3	products PAH sample preparation ples of transfers – TLC separation to preparative column . From TLC to silica gel column From TLC to Sephadex LH-20 Separation systems on two phases –	433 434 435 436 436
	14.2 14.3 Separa 15.1 15.2 15.3 15.4 Enrich 16.1	 14.2 Ethoxylated fatty alcohols

1	React	ion products	445
	1.1	Methyl ricinoleate	445
	1.2	4,4'-Methylene di(phenylisocyanate) diurethane	447
	1.3	2,2,6,6-Tetramethylpiperidine derivative	448
	1.4	Components of two reaction products	449
	1.5	Tetramethylthiuram disulfide mother liquor residue	450
	1.6	Separation of N-vinylformamide and formamide	452
	1.7	Addition of 2,3-epoxypropyl phenyl ether to acetic acid .	453
	1.8	Chlorination of hydroquinone	454
	1.9	Diphenyl ether derivatives	
	1.10	Bis-(4-hydroxyphenyl) sulfone derivatives	458
	1.11	Oligomerization of 4,4'-thiobis-	
		(2-tert-butyl-5-methylphenol)	459
	1.12	Separation of fluorene-bis-9,9-propionitrile (reactant)	
		from fluorene-bis-9,9-propylamine (product)	459
	1.13	Synthesis of 1,4-dibromo-2,3,5,6-tetramethylbenzene	461
	1.14	Synthesis of N-butyl-(2-hydroxyethyl)carbamic ester	462
	1.15	Synthesis of a symmetrical p-xylylene derivative	464

	1.16	New con	npound (Hofmann degradation)	466
	1.17	Unknow		
	1.18	Dipheny	l ether 4,4´-dicarboxylic acid derivative	469
	1.19	Bispheno	ol A derivative	470
	1.20	Synthesis	s of an aromatic phosphorus compound	470
	1.21	Synthesis	s of 2-chloromethylbenzenesulfonamide	471
	1.22	Separatio	on of two by-products of a synthesis	472
	1.23	Synthesis	s of (±)-3-benzyloxy-1,2-propanediol	473
	1.24	Trimethy	vlolpropane acrylates	475
	1.25	Acetylati	on of 2,5-dihydroxybiphenyl	476
	1.26	Addition	of methacrylic glycidyl ester to	
		biphenyl	-4,4´-dicarboxylic acid	477
2	Extrac			479
	2.1	Reaction	product of 2-mercaptobenzothiazole and	
			methane	
	2.2		de determination	
	2.3		liquor residue	
	2.4		monostearate	
	2.5		tion of maleic anhydride	
	2.6	Quantita	tion of b-estradiol in soybean oil	486
	2.7		of a film finishing agent	
	2.8		of synthetic fragrance mixtures	
3			oducts – Methods of purification	
	3.1	Purificat	ion of racemates	491
		3.1.1	(±)-1-Phenyl-1-propanol and	
			(±)-2-phenyl-1-propanol	
		3.1.2	(±)-2-Amino-1-phenylethanol	
		3.1.3	(±)-a-Methylbenzylamine	
		3.1.4	(±)-1-Phenylethanol	
		3.1.5	1-Phenoxy-2-propanol	
		3.1.6	Purity checks on enantiomers	498
	3.2	Isolation	of unavailable substances from mixtures of	
		oligomer		
		3.2.1	Polyethylene glycol oligomers	
		3.2.2	Polytetrahydrofuran oligomers	
		3.2.3	Poly(dimethylsiloxane) components	
		3.2.4	Polyethyleneglycol dimethyl ether oligomers	
		3.2.5	Polyethylene glycol monomethyl ether oligomers	
	3.3		and purification of standards	
		3.3.1	Dalapon	
		3.3.2	Dalapon-methyl	
		3.3.3	Polyethylene glycol alkyl ethers	
		3.3.4	Ethoxyquin	517

		3.3.5	N-Phenylthiourea	518
		3.3.6	Piperonyl butoxide	520
		3.3.7	Purification of 2,6-dimethylphenol	522
	3.4	Monod	isperse polyethylene glycols	523
		3.4.1	Isolation of free polyethylene glycols	
		3.4.2	Purification of polyethylene glycol monomethyl	
			ethers	524
		3.4.3	O-(2-tert-Butyloxycarbonylethyl)dodecylethylene	
			glycol – Purification	
		3.4.4	O,O'-Bis-polyethyleneglycol derivatives	
		3.4.5	Emulsifier	530
	3.5	Mono-	and bifunctional polyethylene glycols (PEG)	532
		3.5.1	O-(2-Aminoethyl)-O´-methyl-polyethylene	
			glycol 750	532
		3.5.2	Poly(ethylene glycol) butyl ether Mn ca. 206	533
		3.5.3	Polyethylene glycol 200 monooctyl ether	533
4	Natu	ral produ	cts	535
	4.1	Curcun	າin	535
	4.2	Oligosa	ccharides	537
	4.3	Conifer	yl alcohol from Siam benzoin	539
	4.4	Essentia	al oils	543
		4.4.1	Aniseed oil (SD)	546
		4.4.2	Siam benzoin resinoid (AE)	549
		4.4.3	Eucalyptus oil (SD)	
		4.4.4	Fennel oil	
		4.4.5	Clove oil – Eugenol	552
		4.4.6	Peppermint oil (SD)	555
		4.4.7	Rose oil	557
		4.4.8	Tea tree oil (bush oil) (SD)	559
		4.4.9	Cedar oil	561
		4.4.10	Analysis of the essential oil of cinnamon bark	564
		4.4.11	Thin-layer chromatograms of various essential oils	566
		4.4.12	Total contents of drugs and preparative separation	n
			of mixtures (extracts)	576
	4.5	Resins -	– Balsams	577
		4.5.1	Siam benzoin	579
		4.5.2	Sumatra benzoin	581
		4.5.3	Peru balsam	
		4.5.4	Tolu balsam – Tolu resinoid	
	4.6	Podoph	ıyllin	
	4.7	(+)-Rut	in	593
5	Plant	extracts		594
	5.1	Rhatan	y root	595

	5.2	Hop cone extracts	. 598
	5.3	Parsley seed extract	. 601
	5.4	Garlic	. 603
	5.5	Tonka extract	. 605
	5.6	Bergamot oil	. 606
	5.7	Paprika powder (condiment)	
	5.8	Black pepper (condiment)	
	5.9	Curry powder (condiment)	
	5.10	Onions	
6	Comn	nercial products	
	6.1	Sorbitan monooleate (Span®80)	
	6.2	Fatty acid ethoxylates	
	6.3 (±)-Glycerol glycide	. 622
	6.4	Pentaerythritol tetraester mixture	
	6.5	<i>a</i> , <i>a</i> ´-Diglycerol	
	6.6	Triglycerol	
	6.7	Polypropylene glycols	. 626
		6.7.1 Dipropylene glycols	. 627
		6.7.2 Tripropylene glycols	. 627
		6.7.3 Tetrapropylene glycol – Di(trimethylolpropane)	. 628
	6.8	Diacetin	
	6.9	Ethoxylated N,N'-bis-(2-hydroxyethyl)-5,5-	
		dimethylhydantoin (mono- and dioleate)	. 631
	6.10	Functionalized oligoethylene glycols	. 633
		6.10.1 12-Amino-4,7,10-trioxadodecanoic acid	
		<i>tert</i> -butyl ester	. 633
		6.10.2 2-[2-(2-Methoxyethoxy)ethoxy]acetic acid	. 634
		6.10.3 3,6,9-Trioxaundecanedioic acid	
7	Deriva	atives of PEG oligomers	. 635
	7.1	Derivatives of PEG oligomers – Ethers	. 635
	7.2	Derivatives of PEG oligomers – Esters	. 644
	7.3	Pentaethylene glycol monodecyl ether	. 648
	7.4	Triethylene glycol monoalkyl ether	. 649
	7.5	Diethylene glycol (DEG) monoalkyl ethers	. 650
	7.6	Ethylene glycol (EG) derivatives	. 650
	7.7	Commercial products	. 650
		7.7.1 Poly-(1,4-butanediol) adipate	. 651
		7.7.2 Poly-(1,4-butanediol)-bis-(4-aminobenzoate)	. 655
		7.7.3 O,O'-Bis-(2-carboxyethyl)dodecaethylene glycol	. 656
8	React	ion products	
	8.1	Cyclohexane-1,4-dicarboxylic acid derivatives	. 658
	8.2	Trifurfuryl tricarballylate	
	8.3	Ethoxylation of bisphenol A	

	8.4	Transesterification oligomerization	. 662
	8.5	Synthesis of cyclic oligomers	. 664
	8.6	Addition products of oxydianiline and glycidol	. 665
9	Trans	fer from a TLC plate to a preparative column	. 667
	9.1	Phenyl glycide ether	. 668
	9.2	2-Benzylaminoethanol	. 669
	9.3	Perfumes – Eaux de toilette	. 670
	9.4	Galingale	. 673
	9.5	Javanese turmeric	. 676
	9.6	Curcuma rhizome	. 678
	9.7	Ginger rhizome	. 678
	9.8	Angelica root	
	9.9	Lovage root	
	9.10	Myrrh – Essential oil and resin	
	9.11	Vanilla extract (vanilla pods)	
	9.12	Oakmoss – Essential oil	
	9.13	Arnica flowers – Arnica tincture	
	9.14	Propolis	
10		anism of separations on Sephadex LH-20	
	10.1	Chemical structure and elution behavior	
	10.2	Test mixtures – Sense and purpose	
	10.3	Separation of vanillin and isovanillin	
11	Plant	constituents	
	11.1	Comparison separations – DCCC and Sephadex LH-20	
		11.1.1 Vegetable acids	
		11.1.2 Glycosides	
		11.1.3 Flavones	
	11.2	Coumarin derivatives	
	11.3	Incense resin – Incense oil	
	11.4	Separation of blackcurrant juice on Sephadex LH-20	
	11.5	Straightforward separations proposed for plant constituents	
12	Cyclo	dextrins	
13		nplishable preparative separation proposals	
	13.1	Purification of reaction products	
	13.2	Purification of 4-isopropylcalix[4]arene	
	13.3	Purification of tridodecylamine	
14	Epoxi	ide additions – Separations	
	14.1	Addition of propylene oxide to <i>tert</i> -butanol	
	14.2	Addition of glycidol to <i>tert</i> -butanol	
	14.3	Addition of glycidol to cetyl alcohol	
	14.4	Addition of glycidol to benzyl alcohol	
	14.5	Addition of glycidol to myristic acid	
	14.6	Addition of glycidol to bisphenol A	

	14.7	Addition of phenyloxirane to methanol	736
	14.8	Addition of propylene oxide to benzyl alcohol	737
	14.9	Addition of ethylene oxide to benzyl alcohol	738
	14.10	Addition of (±)-propylene oxide to methanol	739
15	Isolati	on of phytosterols, tocopherols, and tocotrienols from	
	enrich	ed samples	740
16	Pharm	naceuticals	741
	16.1	Coronary therapeutic agent	742
		Dequalinium disalicylate	
	16.3	O-(2-Hydroxyethyl)rutinosides	744
17		ymerization products	
	17.1	Polyurethane depolymerization	747
	17.2	Reaction product – Polyester depolymerization	747

Appendix

A.1	Reaction products – Separations on Sephadex LH-20	
A.2	Miscellaneous separations from everyday practice	
A.3	3 Test mixtures and separation mechanisms on Sephadex LH-20	
	Molecular structures and elution behavior	
	Concluding comments and outlook	
A.4	Structural formulae and elution behavior of	
	low-molecular weight compounds on a standard	
	5-m Sephadex LH-20/methanol column	
	(\bigstar) , and an alphabetical list of	
	compounds $\underline{www.vogel-buchverlag.de} \rightarrow lnfo Click$	
Remarks		
Glossary		
Dictionary of chromatography		
Company directory		
Bibliography		
Index		

I Introduction

When speaking of preparative column liquid chromatography, it is necessary to distinguish between the separation of synthesis products on a laboratory scale as well as the isolation of pharmaceutically relevant plant constituents on the one hand and separations and purifications of various products on a production scale on the other. Use of the heading 'preparative chromatography' in both cases is inappropriate because the subject matter will only become apparent from the content.

The products of organic syntheses (reaction mixtures and extracts) and plant constituents are generally separated or purified with the aid of low-pressure column liquid chromatography for reasons of cost. Preparative column liquid chromatography becomes an expensive process if the separations or purifications incur high costs compared to classical methods such as distillation, crystallization, or extraction. These higher costs arise from expensive column packing materials and the large quantities of solvent mixtures required. The lack of versatility and low loadability of the separation systems are also cost factors.

Whenever the main emphasis is placed on industrial column liquid chromatography in the relevant monographs, a systematic treatment of this method in the low pressure range often becomes impossible or is perhaps even undesired (?).

These separation methods most commonly used in the gram, or higher gram, ranges have been, and still are, wantonly neglected. The firm persuasion prevailed, and still prevails, that everyone masters this separation technique. This important area of LC is also overlooked because there are no offers of seminars or workshops held by independent experts which would provide attendees with separation proposals.

Another reason is that laboratory staff generally receive a basic training in classical, i.e. manual, column chromatography. Consequently, there is a dearth of qualified experts required for systematic education in preparative LC. The excellent and highly promising techniques of flash chromatography and liquid-liquid partition chromatography alone with a limited number of documented separations will not suffice to solve the problem at hand.

Simple and powerful separation systems of versatile selectivity employing neat solvents as eluents are required for the economical separation of numerous substance mixtures of widely differing composition. The feasibility of achieving this aim is apparent from the many examples presented in this book.

Modern column liquid chromatography is based on the method of adsorption chromatography developed by the Russian botanist Michael Semenovich Tswett and used by him for the separation of vegetable dyes. In 1906 M.S. Tswett reported the separation of chlorophyll with ligroin on a calcium carbonate column. Separation of the components of the green leaf coloring matter with a neat solvent constituted a decisive improvement in the separation procedure. And although Richard Willstätter may have considered the separation of plant dyes by chromatographic adsorption analysis according to Tswett to be unsuitable for preparative work, the method is justifiably regarded as a preparative technique by today's chromatographers.

What do we mean by preparative chromatography?

The term preparative chromatography is applied generally to a separation technique, in liquid (LC) or gas chromatography (GC), which can be used to isolate pure substances from a mixture of natural or synthetic origin in comparatively large, i.e. preparative, quantities (mg to kg).

Why are preparative separations performed?

In general: The amount of substance to be isolated depends upon its intended use.

- □ Structure elucidation requires 30–50 mg of pure substance.
- □ If the substance is to be used as analytical standard then one hundred or more milligrams have to be isolated.
- □ If it is intended to undertake further reactions with the pure substance then isolation has to be performed on a gram scale.
- □ However, if the preparative method is used to obtain pure substances in production quantities, hundreds of grams or kilograms have to be separated or isolated.

Some of the above reasons for performing preparative separations will now be considered in detail.


- Isolation of pure substances for identification purposes, e.g. the main and/or side components of a synthesis or a mixture of natural products. The amount of substance required will also depend upon the equipment available in a laboratory. As a rule, structure elucidation by ¹H-, ¹³C-NMR, FD-MS, IR, and elemental analysis requires 30–50 mg of pure substance.
- 2) Isolation of standard substances, whether main or side products, for analytical methods such as HPLC or GC.
- 3) Isolation of one or more compounds of a reaction mixture (several hundred milligram or gram quantities) for further reactions before performing a synthesis on a large scale.
- 4) Quantification of the individual peaks by gravimetry (a highly accurate method of determination). This provides quantitative information about a mixture of substances even before the qualitative composition of the sample is known. A prerequisite for the success of this method is: one peak = one substance. Precise sample evaluation by mass determination can indicate whether several components are present in roughly equal proportions or whether just one main component (desired or undesired) is present alongside a number of minor components in the reaction mixture. If the compounds are unknown or not available as standards it is not possible to deduce the amount of substance present

from the peak height of peak area in the case of RI or spectrophotometric detection.

- 5) Gentle purification of a reaction product. In this way, time consuming and sometimes difficult purification operations such as distillation, crystallization, re-precipitation, etc. can be replaced by this more effective method. It provides a means of gently removing, e.g., a high-boiling solvent or excess reaction component.
- 6) Straightforward isolation of intermediates in multistep preparations.
- Active principles of plant origin or drugs can readily be quantitatively isolated in large amounts from plant extracts or from pharmaceuticals such as tablets or dragees.
- 8) Isolation or enrichment of side components of a reaction or constituents of natural products present only in small quantities (e.g. <1%).
- 9) Preparative isolation of a substance in gram quantities is often less expensive than optimizing the preparation reaction.

Bearing in mind that it permits separation and isolation at room temperature, the above reasons for carrying out a preparative separation and isolation of pure substances also underline the significance and advantages of this method in comparison with other separation techniques: even sophisticated distillation methods (with scope for variation of temperature, vacuum, high performance columns, etc.) are unable to separate many compounds from their mixture. Nor do other methods of isolating pure substances such as extraction, fractional crystallization, or sublimation always lead to the desired result.

The aims and objectives of preparative LC are summarized in Scheme I.1.

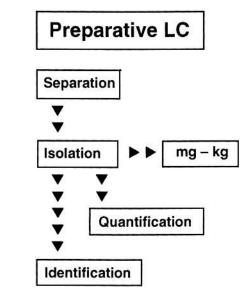

Scheme I.1

The various preparative separation methods in LC

In the course of recent years, various LC methods have been developed for the separation and isolation of pure substances from mixtures of varying degrees of complexity, with the scale-up of analytical HPLC separations to columns of correspondingly greater dimensions (ID, length, coarser-grained stationary phases; 30-60 µm average particle size) probably representing the most common approach. In case of very good resolution, analytical columns can also be overloaded for semipreparative separations. The main drawback of this method is often the insufficient sample capacity (amount of substance/separation = loadability) of the columns in the sense of a truly preparative separation. In addition, insufficient solubility of the sample in the mobile phase can also be an obstacle. It is common knowledge that such adsorption, partitioning, or reversed phase chromatographic separations on silica gel, RP, diol, or CN phases can only be accomplished with organic or organicaqueous solvent mixtures. Consumption of large amounts of solvents, which is unavoidable in preparative column liquid chromatography, means that recovery and regeneration of the binary or even ternary solvent mixtures is far from simple because distillation rarely leads to the desired effect, except in the case of azeotropic mixtures. The "loss" of mobile phase can be limited on use of a peak-controlled fraction collector by collecting that part of the eluate which is free of peaks in a receiver, i.e. that part where the detector signal corresponds to the baseline. In addition to these numerous variants of column chromatographic methods, use is also made of extraction processes (liquid/liquid partitioning) in the form of a chromatography with complex solvent mixtures as stationary and mobile phases. The sample capacity of the last-mentioned methods is generally even lower. The solvent mixtures used are anything but inexpensive in the quantities used and their recovery is far from straightforward. The fraction collectors often have to be cooled to prevent excessive losses of the frequently highly volatile mobile phase. Regardless of whether preparative column liquid chromatography is performed with a solid or liquid stationary phase, each sample mixture will require a different separation system (stationary and/or mobile phase), i.e. within a given separation operation the chromatographic conditions have to be adjusted to the separation task in question, i.e. optimized, because differently composed samples cannot be separated with a single solvent mixture as eluent on one of the above-mentioned stationary phases.

The total absence of doubt concerning the above separation methods is very evident. There are only few users of preparative LC who frequently have to solve very different separation problems and also have access to different separation procedures so that they are able to use that method which holds the greatest promise for a specific problem. There is more demand for straightforward, universally applicable yet powerful separation methods, i.e. there is a wish to be able to solve a variety of separation problems with just one separation system. It is the author's intention to demonstrate, with the aid of numerous examples taken largely from practice, that Sephadex LH-20 with a single or neat volatile solvent as mobile phase largely fulfills these requirements of the practical separation scientist. The following preparative separation methods are currently the ones most commonly used:

- □ Flash chromatography
- □ Low pressure chromatography
- □ Medium pressure chromatography
- □ High pressure chromatography
- Droplet countercurrent chromatography (DCCC)
- □ Rotation locular countercurrent chromatography (RLCC)
- □ Centrifugal countercurrent chromatography (CCC)

Preparative gel chromatographic separations on Sephadex LH-20 with various neat solvents as mobile phases was practiced with great success even several years before introduction of the above preparative methods. The continuing great demand for this separation method is evidenced by the increase in the number of separation units: we operated thirteen complete separation assemblies together with additional sets of columns with various mobile phases and in various dimensions.

Preparative gel chromatography on Sephadex LH-20

Sephadex LH-20 has dual lipophilic (L) and hydrophilic properties (H). This crosslinked dextran gel derives its lipophilic character from the isopropyl groups introduced by hydroxypropylation of Sephadex G-25, which possesses only hydrophilic properties. The hydrophilicity is due to the numerous hydroxyl functions present. As a result of this dual nature, Sephadex LH-20 swells not only in solvents of weak and medium polarity but also in strongly polar ones. The degree of swelling increases with increasing polarity of the solvent (= mobile phase). The principal properties and physical data of the solvents used by the author as mobile phases are listed in Table II.2.

The exclusion limit of Sephadex LH-20 at maximum swelling, such as occurs in water, dimethyl sulfoxide (DMSO), *N*-methyl-2-pyrrolidone, and methanol, is 4000. In other words, all compounds with a molecular mass of \geq 4000 can only pass through the interparticle spaces and thus elute without separation.

It is not so well known that in addition to molecular sieve separations in order of decreasing molecular mass it is also possible to perform highly selective separations with neat solvents as mobile phases. The most effective of these are selective adsorption chromatographic separations from neat solvents of compounds having molecular masses below 1000. Comparable swelling capacities of Sephadex LH-20 in methanol, dichloromethane, and water make the dextran gel a universal column packing material. The range of applications is further extended by the possibility of using other mobile phases differing in their solvent properties and polarities. Depending upon the polarity of the mobile phase and the compounds to be separated it is possible to perform normal phase, reversed phase, and selective adsorption chromatography on hydroxypropylated dextran gel. The most commonly used mobile phases are methanol, acetone, dichloromethane, and water. Less common use is made of chloroform, *N*-methyl-2-pyrrolidone, and ethyl acetate, although the last-mentioned mobile phase forms a highly selective separation system with the dextran gel and is thus ideally suited for exceptional separations.

The most straightforward separation systems are well suited for separating structural and positional isomers, as well as diastereomers and *cis/trans* isomers. The elution or retention times are determined primarily by the heteroelements in the form of functional groups. It would appear very important to emphasize that separations are also possible within the pore volume, i.e. in the "classical" molecular sieve range.

Many years of experience with Sephadex LH-20 have adequately demonstrated that no other separation system displays such a pronounced correlation between chemical and molecular structure and elution or retention behavior as does this straightforward system. This will be considered in greater detail in the individual chapters of this book.

We use the dextran gel in columns of various dimensions. The standard columns or column assemblies have an ID of 25/26 mm and the gel bed height or length ranges from 1 cm to 500 cm. The same gel bed lengths are also used for larger ID columns (52, 63, 70, 75, and 100 mm ID).

The samples to be separated are of various origins, and represent a very wide range of chemical classes. Most common are various reaction products followed by various extracts which are sometimes of highly complex composition.

Successful separations are readily predictable because solubility of the sample in one or more of the above solvents essentially guarantees separation. Another invaluable aid is thin-layer chromatography. Preliminary separations on silica gel with dichloromethane, acetone, or ethyl acetate indicate that excellent separations on a truly preparative scale are possible with the same solvents, even on short gel beds. Neither solvent mixtures nor gradients are required to effect highly selective separations. We shall now consider the advantages of Sephadex LH-20 as column packing material after first having examined the principal arguments against using this dextran gel (see below), which, however, in most cases prove to be without substance.

Why is Sephadex LH-20 so "rarely" used as column packing material in preparative LC?

Although Sephadex LH-20 is used more frequently than is generally recognized as column packing material, its applications are limited mainly to biochemistry and natural products chemistry. The gel is used for special separations or for preliminary or intermediate clean-up steps and with the mobile phases used it does not constitute a universal separation system. However, on use of neat solvents as mobile phases Sephadex LH-20 affords extremely selective separation systems with which we can solve a very wide range of preparative separation problems on a day-to-day basis and which we have used at a constant success level for many years. The main arguments cited against use of Sephadex LH-20 as column packing material will now be briefly presented.

□ Substances can only be separated in order of decreasing molecular size on Sephadex LH-20.

False: With the aid of neat solvents, extremely selective separations of lowmolecular substances from all compound classes can be accomplished which have nothing to do with molecular sieve chromatography (MSC).

- □ Sephadex gels are only suitable as column packing materials for biochemists and natural products chemists.
- □ Handling of soft gels in large diameter and long columns is beset with problems.

False: We have demonstrated the very opposite to be true day by day over several decades.

□ Chromatography in glass columns, especially in large diameter and long columns, is an obsolete method.

False, if the real advantages have been recognized.

Selective separations cannot be performed on short gel bed heights.
 False: We can demonstrate that the opposite is true with the aid of practical examples.

What are the advantages of Sephadex LH-20 over other column packings?

With neat solvents such as methanol, acetone, dichloromethane, chloroform, ethyl acetate, *N*-methyl-2-pyrrolidone as mobile phases the hydroxypropylated dextran gel forms a straightforward, universal, and powerful separating system. If, on the other hand, silica gel or surface modified silica gel is used then selective separation systems will rarely be attained with neat solvents. The three principal characteristic features mentioned so far are recapitulated below:

Straightforward separation system. Only neat, generally volatile solvents are used as mobile phases. Work-up of the eluate fractions in a rotary evaporator (water bath temperature 50–60 °C/ca. 50 mbar) is easy and fast, and above all gentle on the compounds to be isolated.

Universal separation system. It combines several separation mechanisms such as separation according to decreasing molecular size, separation according to differing, sometimes highly selective adsorption effects, while also offering the possibility of partition chromatographic separations. Universal separation systems are particularly important for the practical scientist who very frequently faces different separation problems.

Powerful separation system. These straightforward and universal separation systems are generally highly selective, i.e. they have high loading capacities. Load-abilities of 300 mg or more of sample per gram of dry gel often encountered in practice are indeed very high.

This book has the aim of contributing to our understanding of the separation mechanisms of low-molecular weight organic compounds, and also of inorganic salts, on the hydroxypropylated dextran gel Sephadex LH- 20 with neat solvents as eluents. The numerous examples of separations taken from everyday practice and selected according to various criteria represent just a fraction of the thousands of separations performed. For reasons of patent protection, many extremely interesting separations of reaction mixtures could not be divulged. Instead, separations of comparable or very similar reaction mixtures are shown. As a rule, the test mixtures were selected in such a way that they showed a pronounced similarity to such reaction mixtures. The nature, number, and position of heteroatoms or functional groups responsible for the differences in elution behavior were identical or essentially comparable.

The following colored test chromatogram (see Fig. I.1) gives a foretaste of the separation performance of Sephadex LH-20.

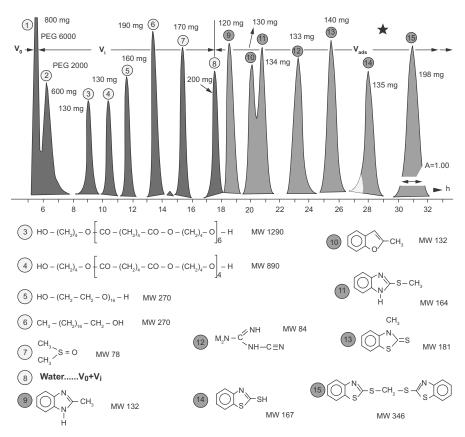


Fig. I.1 Test chromatogram.

A coloured version of this picture you can find in InfoClick on www.vogel-buchverlag.de.

II Theoretical Part

As already mentioned, the term chromatography is used to describe a situation in which the components of a mixture are separated with the aid of a system consisting of two immiscible phases. In the following we shall consider exclusively liquid-chromatographic separation processes performed as column chromatography. Apart from specific chromatographic separation mechanisms such as reversed phase, ion exchange, gel, ion-pair partitioning, and bioaffinity chromatography, chromatographic separations can essentially be described in terms of two physico-chemical processes, viz. **adsorption** and **partitioning**.

A more detailed treatment of adsorption and partition chromatographic separation processes and specific separation processes will be found in relevant monographs.

The chromatographic separation process in a chromatogram is basically described or defined with the aid of the following chromatographic quantities. In column liquid chromatography a chromatogram refers to the elution profile of a substance mixture in the form of peaks recorded on chart paper. The peak shape merely reflects the change in concentration of the eluted substance as a function of time. The individual peaks are described or characterized by the following parameters.

Retention: This includes the retention time (t_R) , the retention volume (V_R) , capacity factor (k), and the relative retention (α) , also known as the separation factor. These parameters are used to describe the peaks of a chromatogram; we use the following schematic representation.

The dead time (t_0) is the time required by a non-retained substance to pass through a column; in other words it is also the time required by the mobile phase. If a compound is retained on passage through the separation system, then it will appear after the elapse of the dead time t_0 , i.e. it shows a **retention time** t_R . This delay or retention comes about because the compounds not only spend time in the mobile phase but also interact with the stationary phase. If two substances are held up to different extents then their elution is staggered (t_{R1} and t_{R2}) and separation occurs as shown Fig. II.1. The **net retention time** t_R' is calculated from $t_R - t_0$.

In order to compare separations in columns of differing dimensions and at different flow rates, the retention time is transformed into a dimensionless quantity, the capacity factor k, in the following way.

$k = \frac{t_{\rm R} - t_0}{t_{\rm R}} = \frac{t'_{\rm R}}{t_{\rm R}}$	(Eq. II.1)
$t_0 \qquad t_0$	

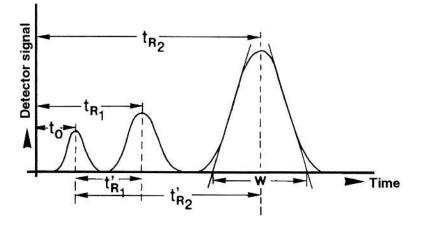


Fig. II.1 Retention times.

k is the ratio of the retention times of a substance in the stationary and in the mobile phase. The equation indicates that the capacity factor k is the net retention time relative to the dead time. Whether a substance resides preferentially in the mobile or in the stationary phase can be described with the aid of the partition coefficient K:

$$K = \frac{Concentration in stationary phase}{Concentration in mobile phase}$$
(Eq. II.2)

The *k* value can also be defined as follows:

$$k' = K \cdot \frac{Volume \text{ of stationary phase}}{Volume \text{ of mobile phase}}$$
(Eq. II.3)

It is thus a logical conclusion that two substances will only be separated if they show a difference.

An appropriate measure is the relative retention α :

$$\alpha = k_2/k_1 \tag{Eq. II.4}$$

Separation is possible only when $k_2 > k_1$. Particularly in the case of preparative separations, the relative retention is the most important measure of the properties of the chromatographic separation system, i.e. of the selectivity.

Before we turn to a further important chromatographic parameter, the resolution *R*, we shall briefly consider a number of further peak characteristics.